

SONICA[®] 2400 S4

REF CODE
Product name
Code
2400 ETH S4 SC V2* 090.027.0003

2400 EP S4 SC V2* 090.027.0004

 SONICA Ultrasonic Cleaner 4.5
 Litres

Accessories

 Stainless steel wire basket
 Internal dim. 265x120x45mm 090.004.0003

Stainless steel lid for 2400 ETH 090.004.0019

 Stainless steel lid with lift sensor
 for 2400 EP 090.037.0002

 Stainless steel auxiliary tank
 SONICA 2400 090.004.0041
 Indicative dim. 300x150x100mm

 Stainless steel beaker holder
 for two glass or plastic beaker 090.004.0036

Plastic beaker 250 ml 090.004.0032

Glass beaker 400 ml 090.004.0033

Further accessories on request

2400 EP S4 SC V2

* *Stainless steel tank made in AISI 304 thickness 08/10 with round tank corners for a better removing of dirt residuals after cleaning process. The tank is conical and realized with a press machine*

General Characteristics

Ultrasonic Cleaning is a safe and efficient modern procedure, which ensures perfect cleaning within the shortest time possible. It has many advantages over ordinary cleaning methods. Ultrasonics can remove even the most tenacious deposits from parts which would be difficult to reach by hand, it also reduces the risk of injury or infection which may occur during manual cleaning.

Notes for use as "Medical Devices"
Denomination: Ultrasonic Cleaning Equipment EMDN: Z12011302

Field of application: Ultrasonic Cleaning Equipment for surgical and dental instruments

Classification: Medical Devices Regulation Class I rule 13, Active MD, Non Invasive MD, Non Implantable MD

SOLTEC S.r.l.
Uffici Commerciali e Sede Legale: Via G. Röntgen 16 - 20136 Milano
 Tel. +39 02 58 308 378 - Fax +39 02 58 308 595 - www.soltec.it - info@soltec.it

Produzione e Magazzino: Via G. Castelbarco, 17 - 20136 Milano - Tel. +39 02 58 324 131

Cod. fisc. e P.IVA: IT11127210158

N. Mecc.: MI 221750 - C.C.I.A.A. Milano N. 1441548

Iscr. Trib. Mi Reg. Soc. N. 343391/8444/41

Cap. Soc. Euro 110.500 i.v. - Iscr. Registro AEE Nr. IT08020000000435



SONICA[®] 2400 S4

Technical specification subject to change

Technical Data

Main line voltage	230/240V~50/60Hz
Main input power/consumption	130 W
Main input power with heating	305 W
Peak HF Ultrasound Power	260 W
US frequency	39 kHz ± 1 kHz with SWEEP SYSTEM Technology
Liquid drain	tube 10x8mm
Transducers	2 (2 Piezoelectric ceramics each)
Tank capacity [L]	4,5
Weight [KG]	3,8
Installation category	Class II according to EN 61010-1
Environmental conditions	Temperature Rel. Hum. 80% up to 31°C with linear decrease up to 50% at 40°C
IP grade	IP 22
CE-compliant	EMC (EN 61326-1) LV (EN 61010-1 EN 61010-2-040) 2017/745/UE Medical Device Regulation Risk Class I according to the rule 13 to the MDR

Overheat switch : Automatic bimetallic Thermostat operating at 140°C fixed on the heating element

Approximately Internal Tank Dimensions


W/D/H (mm)	300 / 150 / 100
Max filling volume tank	4,5 L

External Appliance Dimensions

W/D/H (mm)	325 / 175 / 260
Max operating tank volume	3,15 L

Constructive Materials

Housing	Stainless steel AISI 304
Tank	Stainless steel AISI 304

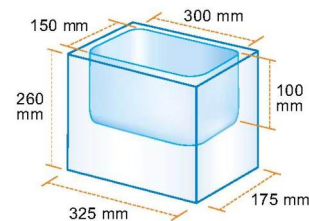

Available Control panel series

ETH

- **Sweep System** generator
- Two line OLED display for quick and intuitive setting of the operating parameters
- Wide setting of time from [1-2-3-4] [5-10-20-25] [10-20-30-40] minutes and cleaning temperature from [35°- 40°] [45°- 50°] [55°- 60°] Centigrades;
- Possibility of remote control functionality activated from control panel (Optional upon request)
- Waterproof control panel



EP

- Sweep System generator
- Two line OLED display for quick and intuitive setting of the operating parameters
- Jet program function to quickly and easily start the cleaning cycle.
- Accelerated heating function and temperature homogenization
- System log function for data traceability
- Pre-set program function
- Lid opening / use control functionality
- Date / Clock function
- Cleaning quality test program function
- Manual free program function
- DEGAS function (0-99 minutes)
- Ultrasonic cleaning function (0-99 minutes)
- Ultrasonic cleaning bath temperature function (up to to 70°C)
- Temperature control via external sensor
- Remote control function (optional on request)



Accessories

External temperature sensor (Optional)



Example of external temperature sensor application



External temperature probe (optional)
product code 091.002.0001

Remote control (optional)



6-pole DIN connector



Multi-core cable with 6-pole (male) DIN connector
(optional) product code 091.002.0002

REMOTE CONTROL technical specifications

In the new SONICA S4 series it is possible to optionally the Remote Control function for both EP and ETH S4 versions. This option can be supplied on request on all models – a feature difficult to find on the market for small and medium-sized “*stand alone*” ultrasonic cleaners.

This function allows to activate/deactivate the ultrasound. Once a temperature has been set on the ultrasonic cleaner it is also possible to turn on and off the heating resistor via remote control.

REMOTE CONTROL mode use and electrical connections

The machines optioned for the Remote Control mode are equipped with a 6-pole (female) DIN connector on the back of the equipment, covered by a protective cap.

A multi-core cable with a 6-pole (male) DIN connector is also supplied with the equipment. The opposite end the cable is unsheathed and individual cables can be connected to the controller.

Control logic and connection to the controller

The main function is the possibility of activating the machine independently of the equipment controller, thus making it possible to decide the start and duration of the cleaning cycle according to the needs.

Switching off the heating resistor is a secondary function, because it always requires the temperature of the cleaning liquid to be set and the heating to be switched on from the equipment panel, as seen in the instructions for use, since the temperature control is always managed by the equipment controller.

This function can be useful if you need to control the temperature of the bath using an external temperature sensor, connected to the controller and immersed in the tank containing the cleaning liquid.

For more detailed information, please refer to the Instruction For Use of our New SONICA Ultrasonic Cleaners S4 series. You will find the connection diagram to the external controller for remote control